Laboratory Washer (Semi-automatic, Compact)

AW47

Capacity

Test tube 450 pcs. (16.5mL)

Easy to use benchtop semi-automatic glass-ware washer.

- Semi-automatic
 - washer, easy to operate by simply setting time then start
- Upward and downward two-way pressurized water jet method with rotating jet nozzles bring high level cleaning. Detergent washing is also available
- With built-in water heater, no boiler piping and water heating system are required
- Optional jet rack is available for hard to clean targets, such as glassware with narrow neck or body

Control Panel






Specifications

Model	AW47		
Cleaning method	Upward and downward two-way pressurized water jet method Rotating jet nozzles (fixed when using jet rack)		
Water supply temp.	Room temp. to 60°C		
Hot water supply system	Heater heating (built-in: 1kW)		
Supply water pressure	0.1 to 0.3MPa		
Glassware stand	Table (Standard), racks (optional)		
Water supply	Open and close by solenoid valve		
Water drain	Natural drainage by water level gap		
Exterior material	Chrome-free electric galvanized steel plate, Chemical-resistant paint		
Interior material	Stainless steel plate		
External dimensions*1	W450×D490×H875mm		
Internal dimensions	W420×D450×H570mm		
Pump	200W		
Turntable	Dia. φ420mm		
Door	Drop down style		
Weight	ht Approx. 43kg		
Power source (50/60Hz)	AC115/220V Single phase with step-down transformer		
	Water supply hose: 2m×1pc. (with coupler), Drain hose: 1.5m×1pc. (I.D.25.4mm)		
Accessories	Vinyl cover 1pc.		
Accessories	Phosphorus-free detergent 1kg (50ml measuring spoon 1pc.)		
	Water faucet		
Consumable	Phosphorus-free detergent		

^{*1} Do not include protrusions.

The length of the power cord is about 2m outside the unit.

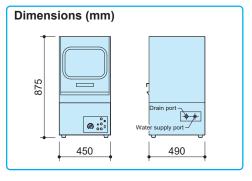

Optional items



Jet rack (glassware not included)



Test tube rack (glassware not included)



Phosphorus-free detergent

No.	Product name	Description	Product code
1	Jet rack	Hold up to 36 pcs. of 100ml flask	291090
2	Test tube rack	Hold up to 450 pcs. of ø18.5mm test tube	291091
3	Detergent	Phosphorus-free detergent 8kg	8190026001